
## Alpinglow Summit ReCover User Manual

The Alpinglow Summit ReCover is a versatile and portable red light therapy system featuring flexible dual-LED panels that deliver 360° full-body treatment. Designed for both relaxation and performance, it can be used in a sleeping-bag style for restorative sessions or opened flat as two separate mats during exercise or stretching.
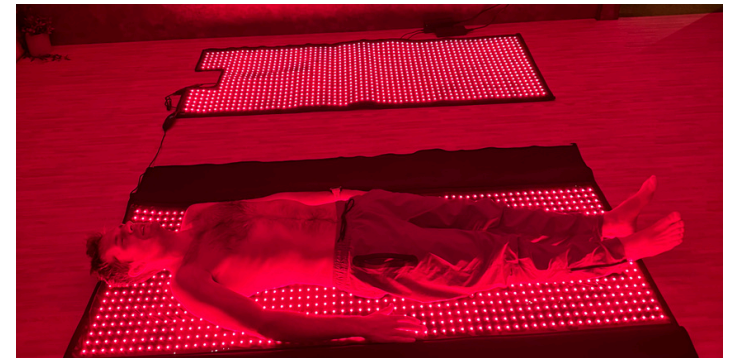
## GETTING STARTED

### Welcome to Alpinglow Summit ReCover – The Future of Light-Powered Recovery and Relaxation

This red light therapy wellness solutions is designed to support recovery, relaxation, and overall well-being.


Designed in Switzerland, the Summit ReCover combines cutting-edge LED light technology with ultimate comfort. It surrounds your body in red and near-infrared light, designed to help you recharge, recover, relax and feel refreshed after daily activity or training.

The ReCover sleeping bag delivers high-quality light exposure through red (660nm) and near-infrared (850nm) wavelengths that are known to support natural processes in the body. This experience can help promote relaxation, improve general vitality, improve sleep and enhance good mood.

The device offers flexible use: you can enjoy it as a full-body sleeping bag for a cocoon-like session, open it fully to create a larger wellness mat, or separate it into two individual mats for versatile light sessions.

For the best experience and results, please read this manual carefully before your first session.

## WHAT'S INCLUDED

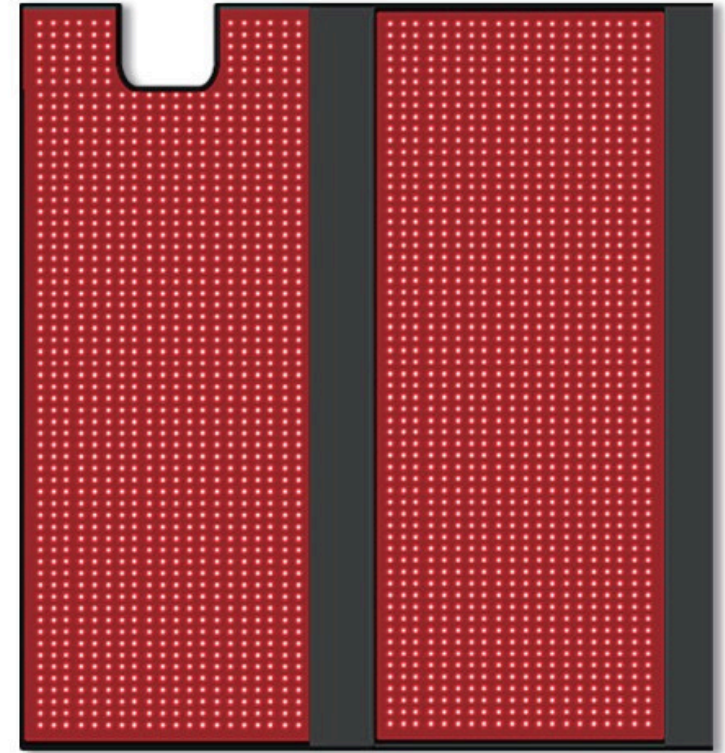1 x Alpinglow Summit ReCover

2 x Controller

2 x Adapter

2 x Power cord

# HOW TO OPERATE

## Step 1

Connect the device to the controller

## Step 2

Connect the controller with the Power adapter, and then plug in the outlet. (Don't use a power bank)

## Step 3

After switching on the mats (remote controls work separately for each part), place the ReCover flat on a dry, smooth, and stable surface. Avoid placing it on very soft or uneven areas, to maintain the longevity of the built-in LEDs.

Choose your preferred mode on each remote control. Lie down and let the light shine directly on your uncovered skin for the most effective experience. Wearing of protective glasses is recommended

Red Light Intensity P1-P5

Timer: 10-90 minutes options

Power button:

1. Power On
2. 10Hz
3. 40Hz
4. Power Off

## Step 4

To set the session time, press the switch button first, then press the timer button on the controller. Each press increases the time in 10-minute increments, up to a maximum of 90 minutes. Press the power button again to confirm and exit the timer setting.

To adjust the light intensity (P1-P5), press the switch button first, then use the + button on the left to increase the intensity or the - button to decrease it. Press the power button again to confirm and stop adjusting the intensity.

\*It is recommended to use it 3-5 times a week, 20-30 minutes each session. Consistent, regular use supports the best long-term wellness and recovery results.

# SMART MULTI-FUNCTION CONTROLLER

Auto shut off Timer. Multiple Mode Options. Pulse mode

Brightness adjustment button  
(P1/P2/P3/P4/P5)



Press and hold the button to turn it ON or OFF

Short press button to switch modes

- 1<sup>st</sup> mode: Constant light mode
- 2<sup>nd</sup> mode: 10Hz Pulse mode
- 3<sup>rd</sup> mode: 40Hz Pulse mode

Time adjustment button

# TECHNICAL INFORMATION & SPECIFICATIONS
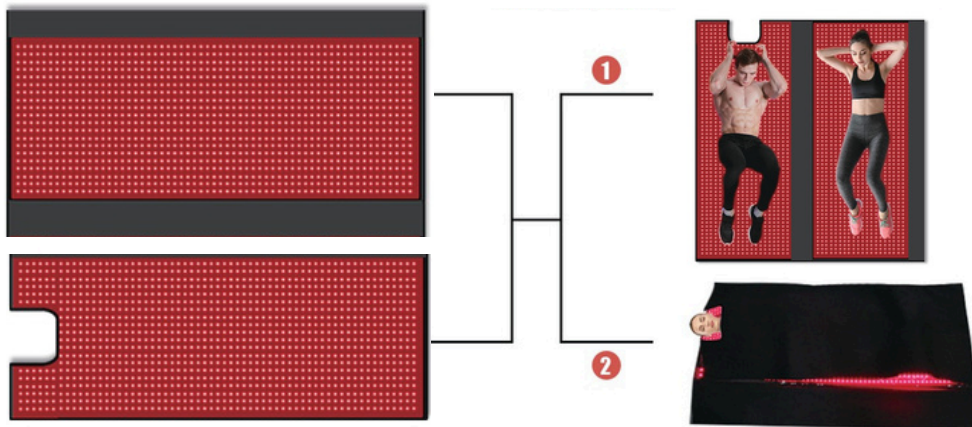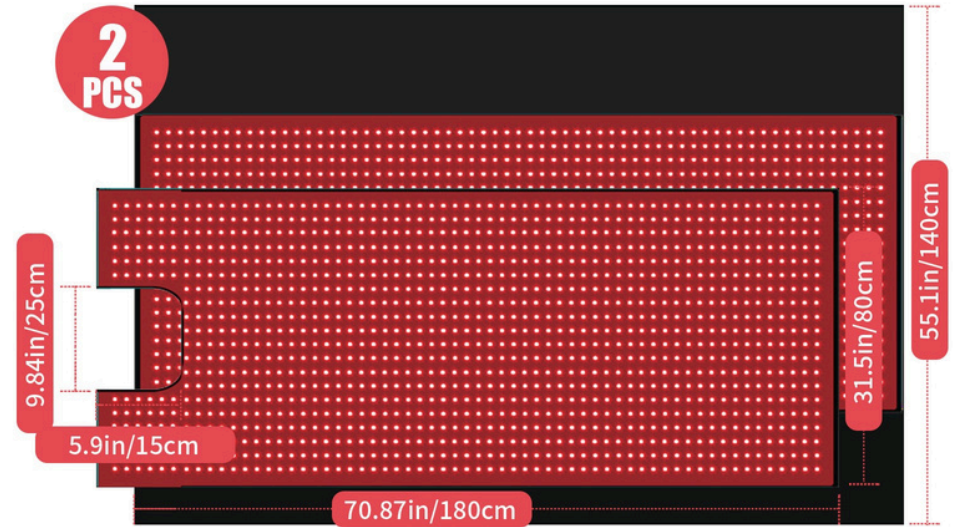
<b>Name</b>	Alpinglow Summit ReCover
<b>Dimensions</b>	180x80cm(sleeping bag) 180x220cm (opened flat in mat style)
<b>Rated watts</b>	270W
<b>Wavelength</b>	660nm/850nm
<b>Number LEDs</b>	2520pcs (660nm:850nm)
<b>Function</b>	Wellbeing

## MORE ABOUT THE DEVICE

Note: The device emits both visible red and invisible infrared light. The infrared light is not visible by human eyes.

5 levels average irradiance (The temperature changes when the intensity is adjusted. The higher the intensity, the higher the temperature.)

Levels \ Average Irradiance	Red Light	Near-Infrared Light (Invisible)
P1	11 mW/cm <sup>2</sup>	15 mW/cm <sup>2</sup>
P2	35 mW/cm <sup>2</sup>	50 mW/cm <sup>2</sup>
P3	65 mW/cm <sup>2</sup>	80 mW/cm <sup>2</sup>
P4	90 mW/cm <sup>2</sup>	111 mW/cm <sup>2</sup>
P5	110 mW/cm <sup>2</sup>	135 mW/cm <sup>2</sup>


Wavelength	Applicable Symptoms
660 nm (Red)	Red light helps relax muscles, supports natural blood circulation, reduces minor aches and swelling, and promotes healthy-looking skin.
850 nm (Near-infrared)	Penetrates deeper into tissues to ease muscle and joint tension, relieve stiffness, and support recovery, relieve Arthritis and tendonitis

## FEATURES:

- High-performance design for full-body wellness sessions, supporting muscle relaxation, comfort, and overall recovery.
- Portable and versatile, allowing use in different settings and situations.
- Equipped with 2,520 LEDs combining 660 nm red light and 850 nm near-infrared light. The near-infrared light is invisible to the eye but can be seen using a phone camera.

3 Chips in 1 Bead
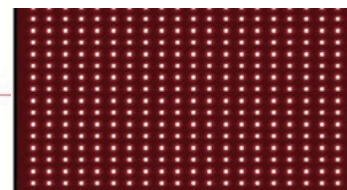
## BENEFITS BACKED BY SCIENCE:

Research shows that the wavelengths in the spectrum of 660-850nm show optimal biotical response.

- Supports weight loss
- Relieves muscle and joint pain
- Reduces wrinkles and signs of aging
- Accelerates wound healing
- Boosts collagen production
- Reduces inflammation
- Promotes hair regrowth
- Decreases pain
- Enhances blood circulation

## HOW TO BEST USE THE SLEEPING BAG:

- For optimal results, it is recommended to use the Summit ReCover 3–4 times per week, with each session lasting around 30 minutes. With consistent use over a few weeks, you can enjoy a noticeable improvement in relaxation, recovery, and overall well-being.
- The ReCover offers three modes: full light normal mode, 10 Hz pulse mode, and 40 Hz pulse mode, along with five energy levels (P1 to P5).
- We generally recommend starting with the full light normal mode at energy level P3, increasing to P5 if you prefer a warmer sensation. For those with heat-sensitive skin or who prefer a gentler experience, the 10 Hz and 40 Hz pulse modes are ideal. In these modes, the light pulses 10 or 40 times per second, creating a rhythmic on-off pattern. This pulsing reduces the perceived heat and provides a more comfortable sensation while still delivering the benefits of red and near-infrared light.
  - 10 Hz pulse mode offers a slower, soothing rhythm that is gentle on the skin and promotes relaxation.
  - 40 Hz pulse mode pulses more quickly, which can enhance circulation and provide a more stimulating experience while remaining comfortable for sensitive skin.
  - These pulse modes allow you to tailor the session to your comfort and preferences without compromising the effectiveness of the therapy.
- TIP Switch on the ReCover sleeping bag 2–3 minutes before use to allow it to warm up slightly, ensuring maximum comfort from the start.

## WARRANTY

Comes with a two-year warranty for reliable use. Kindly follow the instruction in this manual for longevity of the LEDs.

## HOW TO CLEAN:

You must turn off the power first and make sure it is not working then clean it. DO NOT clean and immerse it in water. The red light therapy mat is not waterproof. It only needs to be gently wiped with a dry towel. If there is fog inside the physical therapy mat and cannot see clearly, you can open the zipper and wipe it with a dry towel gently.

## NOTICE:

Although this technology is generally safe and has no known side effects, please follow these precautions:

- This device is not waterproof.
- Do not use the device for more than 30 minutes per session.
- Do not stare directly at the lights.
- The device is not recommended for use while sleeping.

## IMPORTANT NOTICE:

**Do not use the mat on soft or uneven surfaces (such as beds or couches). Always place the device on a firm, flat surface or a stable massage table to ensure proper function and longevity.**

**Avoid excessive bending, folding, or improper handling of the device. Unnecessary or incorrect overbending may damage the internal connections and affect performance.**

## CONTRAINDICATION:

Although red and near-infrared light therapy is generally safe, this product should not be used in the following situations:

- Over open, bleeding, or severe skin injuries.
- By children, pregnant women.
- By people with cancer or suspected cancerous lesions; do not use over malignant or benign tumors.
- By anyone with a history of photosensitive epilepsy, seizures, or loss of consciousness triggered by flashing lights.
- Over recent surgical areas, LASIK or PRK eye surgery, skin filler injections, breast implants, or tattoos in the treatment area.
- By anyone taking medications or using products that may increase light sensitivity, including certain antibiotics, chemotherapy drugs, cold or allergy medications, pain medications, diuretics, citrus essential oils, St. John's Wort, antibacterial soaps, coal tar, or petroleum products.
- By anyone with conditions that may cause sensitivity to light, such as lupus, some forms of eczema, light-induced migraines, albinism, or rosacea.
- By anyone who may be sensitive to red or near-infrared light.

If you have any pre-existing medical conditions or concerns, consult a healthcare professional before use.

## WARNING:

To reduce the risk of burns, fire, electric shock, serious injury, or damage to the product or property, please follow these precautions:

- Only use the provided power adapter and USB cable at the specified voltage. Using other adapters or cables may cause fire or explosion. Do not use the provided adapter or cable to power other devices.
- Do not modify, tamper with, or damage the power adapter, USB cable, or device. Never unplug the power adapter by pulling on the cord.
- This device contains no user-serviceable parts. Do not attempt to disassemble it.
- Inspect the device, remote control, and power cables before each use. Do not use if damaged, malfunctioning, or exposed to water.
- Keep the device, remote, and cables away from flammable or explosive items, chemicals, fire, and heat sources.
- Do not operate the device in the presence of aerosol sprays, corrosive gases, or concentrated oxygen.
- Do not use the device in high temperatures, direct sunlight, under blankets, or under pillows.
- Do not operate while sleeping, driving, smoking, or under the influence of drugs or alcohol. Do not use on sleeping or unconscious people, infants, or children. This device is for adult use only.
- The device is not intended for adults with diminished mental, sensory, or physical capacity unless closely supervised by a responsible adult.
- The 10 Hz pulsed-light mode may affect users with photosensitive epilepsy or other light sensitivities. If you or anyone nearby experiences dizziness, blurred vision, eye or muscle twitching, loss of consciousness, disorientation, involuntary movements, or convulsions, stop using the device immediately and consult a physician.
- Stop using the device immediately if you feel pain or discomfort.

Thank you for investing in your health and well-being with Alpinglow Summit ReCover.

[www.alpinglow.ch](http://www.alpinglow.ch)

[alpinglow\\_switzerland](https://www.instagram.com/alpinglow_switzerland)

[Alpinglow](https://www.facebook.com/Alpinglow)

[info@alpinglow.ch](mailto:info@alpinglow.ch)

+41 79 674 71 07

Crans Montana, Switzerland

Designed and Tested in Switzerland

[www.alpinglow.ch](http://www.alpinglow.ch)

Disclaimer: Red light therapy panels are intended for general wellness purposes only and should not be used as a substitute for professional medical advice or treatment. Individual results may vary. Consult a healthcare professional before use.

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