Alpinglow InfraMerino

Wear Your Wellness.

GETTING STARTED

The Alpinglow InfraMerino is a unique, all-in-one red light therapy beanie with a built-in battery, so you can enjoy its benefits completely cable-free while wearing it. Made from high-quality GOTS-certified merino wool, it is soft, breathable, and naturally odor-resistant. The integrated red and near-infrared light insert delivers wavelengths of 660 nm and 850 nm (invisible light) while remaining completely discreet. This wearable therapy supports healthy hair growth, strengthens hair follicles, improves scalp circulation, and helps reduce inflammation, making it a stylish and practical way to incorporate advanced photobiomodulation into your daily hair care routine.

For best results, use the beanie for 20–30 minutes per session, 3–7 times per week. If you're new to red light therapy, start with 2-minute sessions and gradually increase to 10 minutes over the first two weeks as your body adapts.


Your body benefits from every session. Some users notice changes within days, while improvements in most concerns typically appear after 8–12 weeks of consistent use.

Consistency is essential.

HOW IT WORKS

660 nm red light stimulates cellular metabolism, supports a healthier scalp, and accelerates natural healing processes. This wavelength is known to improve follicle activity and encourage hair regrowth.

850 nm near-infrared light penetrates deeper into the scalp tissues, improving circulation, reducing inflammation, and supporting stronger, healthier hair growth from the roots. Note: 850 nm near-infrared light is invisible to the human eye.

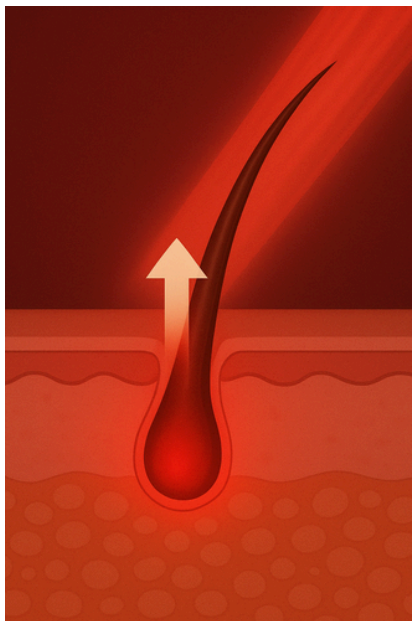Together, 660 nm red light and 850 nm near-infrared light boost follicle activity, improve scalp circulation, reduce inflammation, and create ideal conditions for stronger, healthier hair growth.

OPERATION INSTRUCTION

OFF

No press
->0Hz, Steady Lighting

First short Press ->
Turn ON,
->Level-1 brightness

Frist Press
->Wave Mode,

Second Press,
->Level-2 brightness

Second Press,
->10Hz Pulse Mode

Third Press,
->Level-3 brightness

Third Press,
->40Hz Pulse Mode

Default Timer:
20 Min

Long press -> turn OFF the device.

HOW TO USE

1. Use the charger to fully charge the insert. Do not turn on the device while it is charging. Use a DC 5V 2A power adapter (not included).

2. Put the inner layer in the beanie.
3. Set the Pulse and Brightness Level

Brightness Levels

- Level 1: Low intensity, suitable for relaxation and longer sessions (30 minutes).
- Level 2: Medium intensity, recommended for standard sessions (25–30 minutes).
- Level 3: High intensity, highest irradiance. Recommended for faster results, deeper tissue stimulation, or persistent areas. Use for shorter sessions (10–25 minutes).

Pulse Modes

- Steady Lighting: No pulsing. 0Hz. Provides consistent stimulation, ideal for longer sessions and overall scalp health.
- Wave Mode: Light gradually fades in and out. Provides gentle, comfortable stimulation for the scalp, ideal for sensitive users.
- 10 Hz Pulse: Light flashes 10 times per second. Good for sensitive scalp areas. May help with gentle stimulation of hair follicles.
- 40 Hz Pulse: Light flashes 40 times per second. Supports deeper scalp tissue stimulation and follicle activation. Often linked to neuromodulation benefits.

TECHNICAL PARAMETRES OF RLТ INSERT

- Working Voltage: DC 12V
- Power: 15W
- LED Type: 5050 SMD
- LED Quantity: 150 pcs
- Wavelength Ratio: 660 nm : 850 nm = 1 : 2
- Dimensions: 30 × 19 × 16 cm
- Default Timer: 20 Min
- Dimming Levels: Level 1 to Level 3
- Irradiance: 12mW/cm² (Level 3)

MATERIAL & ORIGIN OF THE MERINO BEANIE

- 100% GOTS-certified organic merino wool (South America)
- Yarn production and dyeing: Schoeller GmbH & CoKG, Hard, Austria
- Knitted in Switzerland (Turgi)

ACCESSORIES

- Merino beanie
- Red light therapy insert with built-in batteries
- USB power cord

MERINO BEANIE

FEATURES

- Breathable, moisture-regulating merino wool
- Naturally odor-resistant
- Durable and comfortable

SUSTAINABILITY

- Compostable natural fiber
- Produced without harmful substances (Detox to Zero)
- Wool from controlled organic animal farming (KBT)
- EXP: chlorine-free and plastic-free finishing

CARE INSTRUCTIONS

- Washing: Wool cycle max. 30°C, low load, spin max. 400 RPM, use wool/delicate detergent only.
- Drying: Lay flat in shade. Not tumble-dryer safe.
- Storage: Store flat; use moth protection if needed.

!!!!Always remove the red light therapy insert before washing the beanie. The insert should never be washed or dried—only the merino beanie itself can be washed.!!!!

RED LIGHT THERAPY INSERT

STORAGE CONDITIONS OF RLT INSERT

- Store the red light therapy beanie and all accessories in a clean, dry location.
- Safe operating temperature: 0°C to 40°C
- Relative Humidity: Below 50%.
- Do not expose the device to freezing temperatures.
- The device is not waterproof. Do not use in rain or snowfall unless fully protected under a waterproof hood or jacket or umbrella.
- Ensure the Alpinglow InfraMerino insert remains dry at all times.

CARE & CLEANING OF RLT INSERT

- Always power off and disconnect the device before cleaning.
- Do not fold or sharply bend the red light insert.
- Clean using a soft, dry or slightly damp microfiber cloth.
- Do not wash or spray water or liquids directly onto the device.

Important:

- Do not use alcohol, chemical cleaners, or abrasive materials.
- Do not operate the device while it is charging

SAFETY INFORMATION

- Do not use this device on infants or young children.
- The beanie is not waterproof. Keep it away from liquids, moisture, and essential oils.
- Only use the power cable for charging provided with the product.
- Do not use this device if you are allergic to neoprene or any of the materials in the beanie.
- Avoid using creams, balms, liniments, or ointments that produce heat before treatment, as this may cause irritation or burns.
- Do not use the beanie on open wounds.
- Avoid use if you are sensitive to light or currently taking medications that increase photosensitivity, as this may cause a reaction. Discontinue use immediately if you feel discomfort or unusual sensitivity.
- Consult a healthcare professional before using the beanie if you are pregnant, have photosensitive conditions, recently received steroid injections, use topical steroids, or have any concerns related to your health or medications.
- Never modify, disassemble, or attempt to repair the beanie.
- Keep the device away from children; it is not a toy.
- Avoid direct eye exposure to the LEDs.
- Do not use the device if it is damaged.

WARRANTY

1. Our products are warranted against defects in materials and workmanship for a period of TWO (2) YEARS from the purchase date.
2. Free return is guaranteed for a period of one MONTH.
3. Warranty does not cover the following damage(s): product been modified, improper installation or use, improper clean, electrical surges, or nature caused damage.
4. Warranty does not cover the loss of third party with direct or indirect damage(s)

DISPOSAL

Do not dispose of this device with household waste. Use a designated collection point for electrical recycling to protect the environment and conserve resources.

DISCLAIMER

This product is not intended to diagnose, treat, or cure any condition. Consult a qualified health professional before use for any medical concerns.

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