

Alpinglow Summit Mat

User Manual

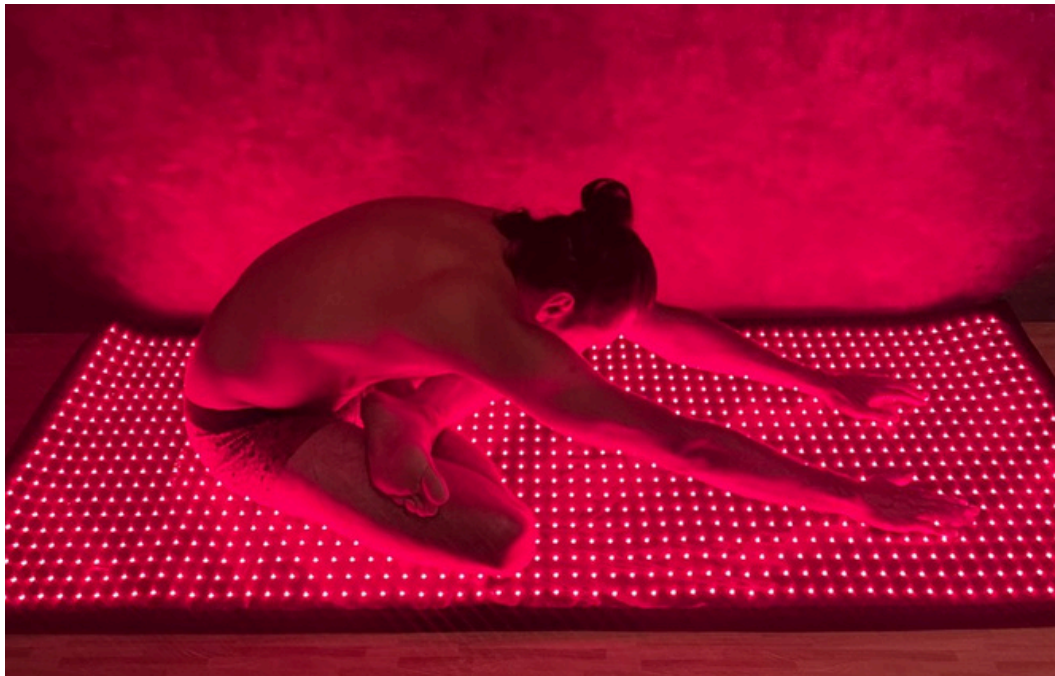The Alpinglow Summit Mat is a versatile, portable red light therapy solution with flexible Dual LEDs, providing full-body therapy and enhancing wellness during exercise or relaxation.

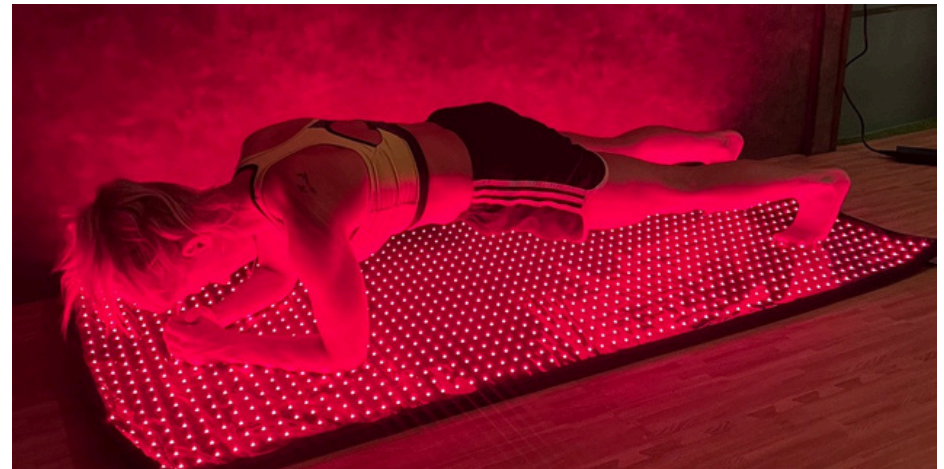Your Red Light Therapy Mat

Before starting your first session, please make sure to read this manual as it may explain important operating procedures. Notes: Each LED inside has three chips . Near Infrared 850nm is invisible to the human eyes.

This package includes

- 1 * Red light Therapy Mat
- 1 * Power Cord
- 1 * Adaptor
- 1 * Connecting Wire

The red light device uses a combo of red 660nm and Infrared 850nm LEDs which is chosen to give the ideal wavelengths found in the scientific research of the last three decades.

660nm red light

This light can be seen by the human eyes and targets skin cells.

850nm near-infrared light


This light cannot be seen by the human eyes and targets deep tissue.

Technical Parametres


- **LED Quantity:** 1280 pcs LEDs
- **Wavelengths:** 660nm : 850nm = 1 : 2 (3-in-1 LED chip)
- **Power Consumption:** 135W
- **Input Voltage:** AC 100–240V, 50/60Hz
- **Timer Setting:** 10–90 minutes adjustable
- **Light Intensity:** 1–5 levels adjustable
- **Pulse Function:** 10Hz & 40Hz
- **Material:** Medical grade TPU film
- **Product Dimensions:** 180 × 80 cm / 70.9" × 31.5"
- **Controller:** 1 controller with digital screen

Red & Near Infrared Light Therapy Mat

Super Large Red Light Therapy Mat

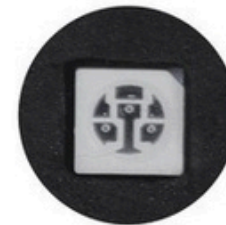3 Chips in 1 Bead

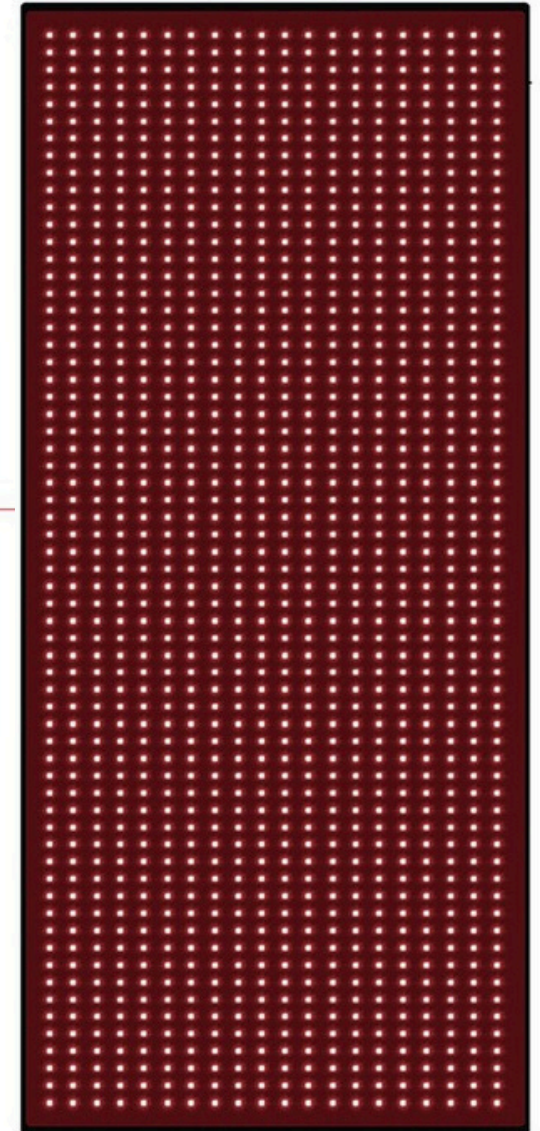How to Operate

Step 1

Connect the device to the controller

Step 2

Connect the controller with the Power adapter, and then plug in the outlet. (Don't use a power bank supply to charge the device. Cause it wildecrease the power and the intensity of the device)

Step 3

After power on, the mat should flat on the dry and smooth floor, then user need wear the anit-red light glasses, choose the mode you like and enjoy the therapy time.

Red Light Intensity P1-P5

Timer

10-90 minutes options

Power

Power On, 10Hz,40Hz,Power Off

Step 4

To choose different timer for each session, press the switch button first then press the timer button on the controller. This will add time in 10-minute increments up to 90 minutes. And press the power button again to turn off the timer function. To choose different intensities, press the switch button first then press the + button for higher intensity, or press the - button for lower intensity. And press the power button again to stop adjusting the intensities. (It is recommended to contact with your skin directly)

How to Best Use the Mat:

It is recommended to use the red light therapy mat minimum 3–4 times per week, with each session lasting 20-30 minutes. After a few weeks of consistent use, you should notice positive results. The mat offers three modes: full light normal mode, 10Hz pulse mode, and 40Hz pulse mode, as well as five energy levels (P1 to P5). We generally recommend starting with the full light normal mode and setting the energy level to P3 - P5 if more heat is desired. For those with heat-sensitive skin, the 10Hz and 40Hz pulse modes are ideal. In these modes, the light pulses 10 or 40 times per second, creating a flashing effect that provides a gentle buffer for sensitive skin.

The 10 Hz pulse mode provides a slower, soothing rhythm that is gentle on the skin and promotes relaxation. The 40 Hz pulse mode offers a faster pulse, which may help stimulate circulation while remaining comfortable for sensitive users. If heat sensitivity is a concern, the 40Hz mode is particularly suitable.

TIP: Switch on the mat 2–3 minutes before use to allow it to warm up slightly, ensuring maximum comfort from the start.

How to Clean:

You must turn off the power first and make sure it is not working then clean it. DO NOT clean and immerse it in water. The red light therapy mat is not waterproof.

It only needs to be gently wiped with a dry towel. If there is fog inside the physical therapy mat and cannot see clearly, you can open the zipper and wipe it with a dry towel gently.

IMPORTANT NOTICE:

Do not use the mat on soft or uneven surfaces (such as beds or couches). Always place the device on a firm, flat surface or a stable massage table to ensure proper function and longevity.

Avoid excessive bending, folding, or improper handling of the device. Unnecessary or incorrect overbending may damage the internal connections and affect performance.

Benefits:

Research suggests that wavelengths in the 660–850 nm range are associated with optimal biological response and are widely used in wellness and recovery applications.

Regular use of red and near-infrared light may support:

- General well-being and relaxation
- Muscle relaxation and comfort after physical activity
- Temporary relief of muscle and joint discomfort
- Improved appearance of skin and support of collagen production
- Reduction in the appearance of fine lines and wrinkles
- Support for normal blood circulation
- Support for natural recovery processes of the body
- A more refreshed and revitalized feeling
- Support for healthy-looking hair
- Reduction of everyday tension and stiffness
- Enhanced post-workout recovery and performance support

CONTRAINDICATION:

Although red and near-infrared light therapy is generally safe, this product should not be used in the following situations:

- Over open, bleeding, or severe skin injuries.
- By children, pregnant women.
- By people with cancer or suspected cancerous lesions; do not use over malignant or benign tumors.
- By anyone with a history of photosensitive epilepsy, seizures, or loss of consciousness triggered by flashing lights.
- Over recent surgical areas, LASIK or PRK eye surgery, skin filler injections, breast implants, or tattoos in the treatment area.
- By anyone taking medications or using products that may increase light sensitivity, including certain antibiotics, chemotherapy drugs, cold or allergy medications, pain medications, diuretics, citrus essential oils, St. John's Wort, antibacterial soaps, coal tar, or petroleum products.
- By anyone with conditions that may cause sensitivity to light, such as lupus, some forms of eczema, light-induced migraines, albinism, or rosacea.
- By anyone who may be sensitive to red or near-infrared light.

If you have any pre-existing medical conditions or concerns, consult a healthcare professional before use.

WARNING:

Safety Precautions

Although this technology is generally well tolerated, please follow these safety guidelines at all times:

- Use only the provided power adapter and USB cable at the specified voltage. Using incompatible accessories may result in damage, fire, or explosion.
- Do not use the provided adapter or cable to power other devices.
- Do not modify, tamper with, or damage the power adapter, cable, or device.
- Always unplug the device by holding the plug—never pull on the cable.
- This device is not waterproof. Keep it away from water and moisture.
- Inspect the device, remote control, and cables before each use. Do not operate if any component is damaged, malfunctioning, or has been exposed to liquids.
- This device contains no user-serviceable parts. Do not attempt to open or disassemble it.
- Keep the device away from flammable or explosive materials, chemicals, fire, and heat sources.
- Do not use the device in the presence of aerosol sprays, corrosive gases, or concentrated oxygen.
- Avoid operating the device in high-temperature environments, direct sunlight, or when covered (e.g. under blankets or pillows).
- Do not look directly into the light source. Always use appropriate protective goggles.
- The device is not intended for use by individuals with reduced mental, sensory, or physical capabilities unless supervised by a responsible adult.
- The 10 Hz pulsed light mode may affect individuals with photosensitive epilepsy or light sensitivities.
- If you or anyone nearby experiences dizziness, blurred vision, eye or muscle twitching, disorientation, involuntary movements, loss of consciousness, or convulsions, discontinue use immediately and consult a physician.
- Stop using the device immediately if you experience pain or discomfort.

Thank you for investing in your health and well-being with Alpinglow Summit Mat.

www.alpinglow.ch

[alpinglow_switzerland](https://www.instagram.com/alpinglow_switzerland)

[Alpinglow](https://www.facebook.com/Alpinglow)

info@alpinglow.ch +41 79

674 71 07 Crans Montana,

Switzerland

Designed and Tested in Switzerland

www.alpinglow.ch

Disclaimer: Red light therapy panels are intended for general wellness purposes only and should not be used as a substitute for professional medical advice or treatment. Individual results may vary. Consult a healthcare professional before use.

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